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Gamma
TECHNICS FOR HEALTH

ELECTRONIC BABY SCALES INSTRUCTION BOOK

Thank you for choosing our Electronic Baby Scales Gamma®. They will help you monitor the growth of your baby, accurately and reliably quantifying how much milk has been assimilated at each feed, taking up to 20 kg, and with 10g divisions. The TARA function allows you to weigh the baby using a soft cloth or towel. The scales come with a useful WEIGHTBLOCK function that allows a steady weight to be displayed despite the child wriggling about.

We are sure that having appraised worthily the high quality and reliability of this device you will become a regular user of the products of Gamma Trade Mark.

Before starting to use the scales please study the instruction book carefully. Following the rules stated in the instruction book will provide for the faultless operation of the scales for a long period of time.

ATTENTION!

- Handle your scales with care as a precision instrument.
- Avoid to expose them to extreme temperatures, damp, knocks, dust and direct sunlight.
- Keep the scales away from sources of heat.
- Take care that no liquids penetrate into the scales, as they are not waterproof.
- Do not leave or use the scales in very damp surroundings, humidity above 85% or contact with water and other liquid parts could cause corrosion and jeopardize the precision of the scales.
- Do not drop the scales or knock on them.
- Clean the platform with a damp cloth. Clean the bottom of the scales with just a soft dry cloth. To clean the scales never use chemical or abrasive products.


INSTALLING BATTERIES

Use a 9V alkaline battery only. Remove the battery cover on the underside of the scales. Insert the battery and refit the cover. Do not throw the batteries away together with your normal domestic rubbish. Take them to a used battery bin.

WEIGHING

Place the scales on a flat and hard surface. A soft and uneven floor(e.g. carpets, linoleum) could cause an alteration to the weight reading.

FUNCTION KEYS

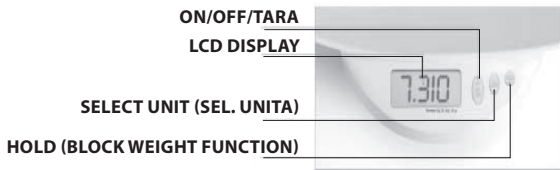ON/OFF/TARA

The ON/OFF/TARA button is for switching the scale on and off and for using the TARA function. Use the TARA function if you intend to place the baby on a towel or soft cloth. In this case, switch on the scales by pressing the ON/OFF/TARA button. The display shows "8888". Wait until the display shows "0.00" then place the towel or cover on the plate. Press the ON/OFF/TARA key and wait until the display shows "0.00"; the scales are now ready for use. To switch off the scales press the ON/OFF/TARA button, otherwise the scales will automatically shut off after 30 seconds.

SELECT UNIT (SEL. UNITA)

Press this key to select the unit of measure you require: kg or pounds. Those scales have a capacity of 20 kg, division 10 g or 44.40 pounds, division 1/2 oz. After switching on the scales, if the display shows 0.00 and a triangle at the top left of the display the weighing will be in kilos and grams; if the display shows 0:0 and a triangle at the bottom left of the display the weighing will be in pounds and ounces. To select the required unit of measurement, switch on the scales and press the SEL.UNITA button.

HOLD (BLOCK WEIGHT FUNCTION)

Due to the baby moving, the weight shown on the display changes all the time. When the baby is laid on the plate, press the "HOLD" key. The display shows "hold" and the scales automatically cal-

culate the average of the weights they have taken. The display then shows the baby's weight and at the top the symbol "△" starts flashing.

SYMBOLS

- "LO" indicates that the batteries are flat. Replace them.
- "ERR" shows the weight is too much.
- The scales are ready for use when the display shows "0.00"; if the instructions have not been followed correctly, the display shows "ERR".

TECHNICAL FEATURES

- Range of measurement: 10 g - 20 kg
- Division: 10 g
- Battery: 1 alkaline battery 6LR61 series
- Working conditions: 0° + 50°C, RH85%
- Storage conditions: 0° + 50° RH 85%
- Dimensions: 32 x 52 x 8 cm
- Weight: 2.20 kg
- Accuracy and repeatability with tolerance level within: 3%

ATTENTION!

Read carefully the instructions for use. Technical alterations are possible.

WARRANTY

Your Electric Baby Scales Model MD6141 are warranted for 2 years from date of purchase. This warranty does not include the batteries.

The warranty does not apply to damage caused by improper handling, accidents, not following the operating instructions or self-maintained alterations made to the device.

The manufacturer shall not be liable for any damages to people, property or pets, either directly or indirectly, caused by the noncompliance with the instructions laid down in the instruction manual.

The warrantee is only valid upon presentation of the warrantee card which was correctly filled in and sealed.

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